Relaxknow

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ENGLISH

| | |
|---------------------------|----|
| Technical Data | 10 |
| Structure and Accessories | 10 |
| Warnings and Precautions | 11 |
| Fault Diagnosis | 11 |
| Usage | 12 |
| Notices / Infos | 14 |
| Clean & Maintain | 14 |
| Disposal | 15 |

TECHNICAL DATA

Power: 25 W

Dimensions: L 390 mm × W 340 mm × H 185 mm

Net weight: 2,4 Kg


Voltage: AC100-240 V 50 / 60 Hz - DC12 V, 2.5 A

Scope of delivery: Massage device, AC adapter, car charger, manual

Standard working time: 10 minutes

STRUCTURE AND ACCESSORIES

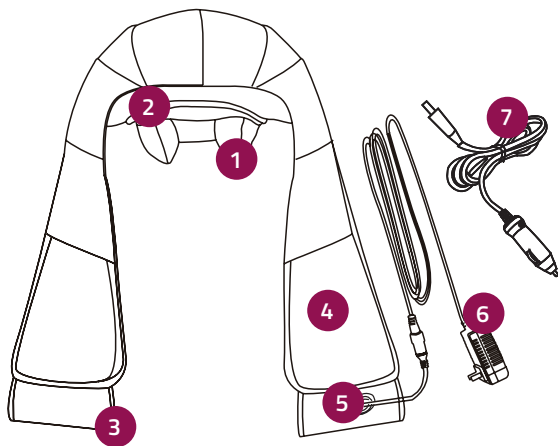1- Removable cover (also washable) with zipper. Pleasant net fabric.

2- 8 massage balls

3- Strap

4- PU Leather

5- DC Connectors

6- Adapter

7- Car Charger

WARNINGS AND PRECAUTIONS

Please read these operating instructions carefully before using the device and keep them in a safe place. Use the device only as described and familiarize yourself with all operating and safety instructions before using the product.

- The following people should not use this device: People with osteoporosis, heart pacemakers, heart disease, pregnant women, people with wound injuries and children under 12 years of age.
- Do not use the unit if it is wet or damp.
- Avoid pulling the power supply unit.
- Dismantling or rebuilding the device can cause damage and injury.
- Use the device for approx. 10 minutes and then leave it switched off for approx. 20 minutes.
- If the power supply or adapter is damaged, please stop using it immediately.

FAULT DIAGNOSIS

- If the device cannot be switched on, please check that the plug is correctly connected and check the power supply.
- If the noise of the device becomes louder and louder, it may be because it has been in continuous use for too long. Turn off the unit for at least 30 minutes before restarting.
- If the unit turns off during use, it may be because the engine is overheated. Leave the unit switched off for at least 30 minutes before restarting.

The device is recommended for the following people:


- People who sit/stand in the same position for a longer time (for example, sitting at a desk in the office for a long time)
- Nervous, tense people
- People with insomnia
- People with stiff neck or stiff shoulder area
- People with headaches or dizziness

You can use the Massage Device for the whole body as shown in the picture:

USAGE

Display panel overview:

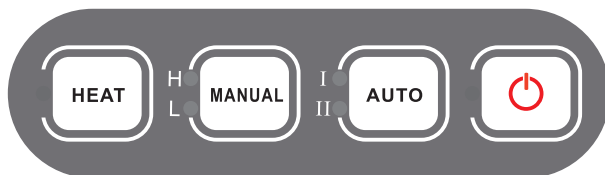How to use the Massage Device:

- Plug the AC power cord connector into the main interface of the unit (make sure it is plugged in to prevent it from falling down).
- Plug the AC adapter into the wall outlet. If the indicator light flashes red, the unit is correctly supplied with power.
- Place the Massage Device on the body area of your choice
- Place your hands in the arm straps and keep them taut when applying to the neck, shoulder, back or hip.

SWITCH ON MASSAGE DEVICE Press the On/Off button on the control panel. The red control lamp lights up.

MANUAL MODE Manual button to adjust the speed and intensity. For a higher speed the „H“ for „Highspeed“ flashes red. For less speed the „L“ for „Lowspeed“ lights up.

AUTO MODE Press the „Auto“ button to enter automatic mode. The unit switches off automatically after 10 minutes.

HEAT MODE The massage device is equipped with a warming function. If necessary, proceed as follows: Press the warming function button.

SWITCH OFF THE MASSAGE DEVICE Press the On/Off button on the control panel.

NOTICES / INFOS

- The device is for private use only, not for medical or commercial use.
- It is suitable for the whole body massage:
- Neck, shoulder, back, legs, stomach, feet and arms.
- If you have medical concerns, talk to your doctor or health care professional before using the Massage device.
- The massage should always be perceived as pleasant. If you find it uncomfortable or even painful, please change the intensity of the massage and remove some tension from the arm straps.
- Do not use the Massage Device while driving.
- Thanks to the car adapter, you can also use the Massage Device on a long trip as a co-driver or at the rest area, for example.
- After 10 minutes of uninterrupted massage, the device switches off automatically.

REINIGUNG & PFLEGE

Attention: Before cleaning the device, make sure that it is switched off and that the AC adapter is disconnected from the wall socket!

This has to be kept in mind:

- Always use a dry cloth for cleaning.
- Never allow any liquid to penetrate into the unit or accessories as this could cause a short circuit.
- Smaller stains should only be removed with a cloth or damp sponge (with a little liquid mild detergent if necessary).
- Cleaners containing solvents should not be used.
- Do not use the Massage Device until it is completely dry.
- You yourself may only carry out minor cleaning work on the device. Do not repair the unit yourself in case of malfunction. Otherwise, any warranty claim would expire.

Due to the practical zipper, the cover can be easily removed and washed at a temperature of 30°C in the washing machine.

The washable cover is inappropriate for ironing.

The washable cover is inappropriate for bleaching.

DISPOSAL

Please dispose the package accordingly. Add carton and cardboard to the waste paper and foils in the recyclables collection.

Applicable in the EU and other European countries with systems for separate collection of recyclable materials.

Do not dispose used appliances in household waste. Every consumer is legally obliged to hand in old appliances separately from household waste (e. g. at a collection point). This is the only way to ensure that old appliances are recycled properly and do not have a negative impact on the environment.

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We, InnovaMaxx GmbH, herewith declare that the model: RXM500 complies with the appropriate basic safety and health requirements of the EU Directive based on its design and type, as brought into circulation by us.

The EMC Directive 2014/30/EU

EN 55014-1:2006+A2: 2011, EN 55014-2: 2015
EN 61000-3-2: 2014, EN 61000-3-3: 2013

2011/65/EU Restriction of Hazardous Substance

IEC 62321-3-1: 2013 IEC 62321-5: 2013 IEC 62321-4: 2013
IEC 62321-7-1: 2015 IEC 62321: 2008 ANNEX C

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For the compilation of the technical documentation is authorized: InnovaMaxx GmbH.

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Innovamaxx GmbH